

Anti Flag-Tag-HRP Mouse Monoclonal Antibody

Catalog Number	Vial Size
LK8002	100 µl
LK8002L	1 ml



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
标准 高效 稳定 Precision Efficient Stable

Market	400-621-0003 marketing@sungenebiotech.com
Support	022-66211636-8024 techsupport@sungenebiotech.com
Web	www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Applications
2E4	Mouse IgG	WB

Description

The DYKDDDDK peptide is a small component of an epitope which does not appear to interfere with the bioactivity or the biodistribution of the recombinant protein. It has been used extensively as a general epitope tag in expression vectors. As a member of Tag antibodies, Flag-tag antibody is the best quality antibody against DYKDDDDK in the research. Flag-tag antibody can recognize C-terminus, internal, and N-terminus Flag-tagged proteins.

Immunogen: A synthetic peptide DYKDDDDK coupled to KLH.

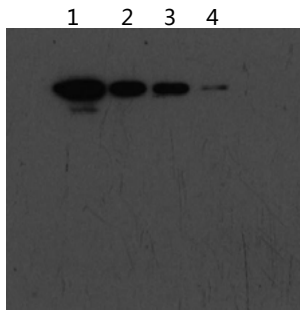

Western blot analysis of Flag-tag fusion protein with anti Flag-tag-HRP mouse mAb (2E4) diluted at 1: 1000. The loading quantity of Flag-tag fusion protein is 100ng(Lane1), 50ng(Lane2), 25ng(Lane3), 12.5ng(Lane4).

Product Information

Recommended dilutions: WB: 1:1000~1:10000,
Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid

Concentration: 1 mg/ml in PBS, pH 7.4, containing 0.02% sodium azide and 50% glycerol

Storage: Store at -20°C or below. Avoid repeated freezing and thawing cycles.

Related Products:

Cat.No.	Name
KM8001	Anti His-Tag Mouse mAb
KM8002	Anti Flag-Tag Mouse mAb
KM8003	Anti c-Myc-Tag Mouse mAb
KM8004	Anti HA-Tag Mouse mAb
KM8005	Anti GST-Tag Mouse mAb
KM8006	Anti V5-Tag Mouse mAb
KM8007	Anti VSV-G-Tag Mouse mAb
KM8008	Anti MBP-Tag Mouse mAb
KM8009	Anti GFP-Tag Mouse mAb
KM8010	Anti T7-Tag Mouse mAb
KM8011	Anti CBP-Tag Mouse mAb
KM8012	Anti HSV-Tag Mouse mAb
KM8013	Anti E2-Tag Mouse mAb
KM8014	Anti S1-Tag Mouse mAb
KM8016	Anti Tag-100-Tag Mouse mAb
KM8017	Anti mCherry-Tag Mouse mAb
KM8018	Anti HPC4-Tag Mouse mAb
KM8019	Anti SUMO-Tag Mouse mAb
KM8020	Anti RFP-Tag Mouse mAb
LK8001	Anti His-Tag-HRP Mouse mAb
LK8003	Anti c-Myc-Tag-HRP Mouse mAb
LK8009	Anti GFP-Tag-HRP Mouse mAb
LK2003	Goat anti-Mouse IgG(H+L)-HRP
LK2004	Goat anti-Mouse IgG(H+L)-AP
GM200G-02	FITC Goat anti-Mouse IgG(H+L)
GM200G-09	PE Goat anti-Mouse IgG(H+L)
GM200G-37	IF555 Goat anti-Mouse IgG(H+L)
GM200G-38	IF488 Goat anti-Mouse IgG(H+L)
GM200G-40	IF647 Goat anti-Mouse IgG(H+L)
GM200G-43	IF594 Goat anti-Mouse IgG(H+L)

For Research Use Only.