

## Anti c-Myc-Tag Mouse Monoclonal Antibody

Catalog Number	Vial Size
KM8003	100 µl
KM8003L	1 ml



天津三箭生物技术股份有限公司  
Tianjin Sungene Biotech Co., Ltd.  
精准 高效 稳定 Precision Efficient Stable

**Market** | 400-621-0003  
marketing@sungenebiotech.com

**Support** | 022-66211636-8024  
techsupport@sungenebiotech.com

**Web** | www.sungenebiotech.com

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.  
This product is guaranteed up to one year from purchase.


### Purified Antibody Characterization

Clone	Isotype	Applications
22E8	Mouse IgG	WB, IF


### Description

c-Myc-tag antibody is part of the tag series of antibodies, the best quality in the research. The immunogen of c-Myc tag antibody is a synthetic peptide corresponding to residues 410-419 of the human p62 c-Myc protein conjugated to KLH. c-Myc- tag antibody is suitable for detecting the expression level of c-Myc or its fusion proteins where the c-Myc tag is terminal or internal.

**Immunogen:** A synthetic peptide EQKLISEEDL coupled to KLH.

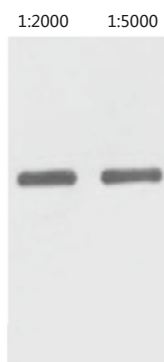

Western blot analysis of c-Myc-tag fusion protein with anti c-Myc-tag mouse mAb (22E8) diluted at 1: 2000(left) and 1: 5000(right).

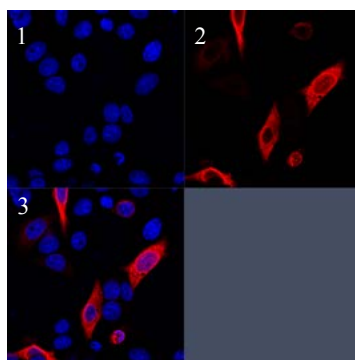

Immunofluorescent analysis of the 293T cell transfected with anti c-Myc-tag mouse mAb (22E8)  
Fig.1: Nuclei stained with DAPI  
Fig.2: Immunofluorescent analysis with c-Myc-tag  
Fig.3: Merge.

### Product Information

**Recommended dilutions:** WB: 1:1000~1:10000,  
IF: 1:200~1:1000

Optimal dilutions/concentrations should be determined by the end user.

**Form:** Liquid

**Concentration:** 1 mg/ml in PBS, pH 7.4, containing 0.02% sodium azide and 50% glycerol.

**Storage:** Store at -20°C or below. Avoid repeated freezing and thawing cycles.

**Related Products:**

Cat.No.	Name
KM8001	Anti His-Tag Mouse mAb
KM8002	Anti Flag-Tag Mouse mAb
KM8004	Anti HA-Tag Mouse mAb
KM8005	Anti GST-Tag Mouse mAb
KM8006	Anti V5-Tag Mouse mAb
KM8007	Anti VSV-G-Tag Mouse mAb
KM8008	Anti MBP-Tag Mouse mAb
KM8009	Anti GFP-Tag Mouse mAb
KM8010	Anti T7-Tag Mouse mAb
KM8011	Anti CBP-Tag Mouse mAb
KM8012	Anti HSV-Tag Mouse mAb
KM8013	Anti E2-Tag Mouse mAb
KM8014	Anti S1-Tag Mouse mAb
KM8016	Anti Tag-100-Tag Mouse mAb
KM8017	Anti mCherry-Tag Mouse mAb
KM8018	Anti HPC4-Tag Mouse mAb
KM8019	Anti SUMO-Tag Mouse mAb
KM8020	Anti RFP-Tag Mouse mAb
LK8001	Anti His-Tag-HRP Mouse mAb
LK8002	Anti Flag-Tag-HRP Mouse mAb
LK8003	Anti c-Myc-Tag-HRP Mouse mAb
LK8009	Anti GFP-Tag-HRP Mouse mAb
LK2003	Goat anti-Mouse IgG(H+L)-HRP
LK2004	Goat anti-Mouse IgG(H+L)-AP
GM200G-02	FITC Goat anti-Mouse IgG(H+L)
GM200G-09	PE Goat anti-Mouse IgG(H+L)
GM200G-37	IF555 Goat anti-Mouse IgG(H+L)
GM200G-38	IF488 Goat anti-Mouse IgG(H+L)
GM200G-40	IF647 Goat anti-Mouse IgG(H+L)
GM200G-43	IF594 Goat anti-Mouse IgG(H+L)

**For Research Use Only.**