

Anti STREP-Tag Mouse Monoclonal Antibody

Catalog Number

KM8015



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
标准 高效 稳定 Precision Efficient Stable

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Applications
6C2	Mouse IgG	WB

Description

Monoclonal antibodies specific to STREP-tag can greatly improve the effectiveness of several different kinds of immunoassays, helping researchers identify, detect, and purify STREP-tag fusion proteins in bacteria, insect cells, and mammalian cells. Anti STREP-tag mouse mAb recognizes STREP-tag placed at N-terminus, C-terminus of fusion proteins.

Immunogen: A synthetic peptide WSH PQFEK coupled to KLH.

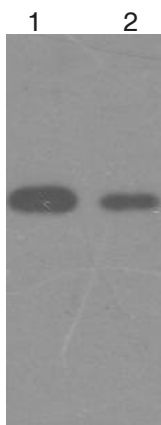

Western blot analysis of STREP-tag fusion protein with anti STREP-tag mouse mAb (6C2) diluted at 1: 1000.
lane1: Protein at 6 µg.
lane2: Protein at 3 µg.

Product Information

Recommended dilutions: WB: 1:1000 ~1:10000,
Optimal dilutions/concentrations should be determined by the end user.

Form: Liquid

Concentration: 1 mg/ml in PBS, pH 7.4, containing 0.02% sodium azide and 50% glycerol

Storage: Store at -20°C or below. Avoid repeated freezing and thawing cycles.

Related Products:

Cat.No.	Name
KM8001	Anti His-Tag Mouse mAb
KM8002	Anti Flag-Tag Mouse mAb
KM8003	Anti c-Myc-Tag Mouse mAb
KM8004	Anti HA-Tag Mouse mAb
KM8005	Anti GST-Tag Mouse mAb
KM8006	Anti V5-Tag Mouse mAb
KM8007	Anti VSV-G-Tag Mouse mAb
KM8008	Anti MBP-Tag Mouse mAb
KM8009	Anti GFP-Tag Mouse mAb
KM8010	Anti T7-Tag Mouse mAb
KM8011	Anti CBP-Tag Mouse mAb
KM8012	Anti HSV-Tag Mouse mAb
KM8013	Anti E2-Tag Mouse mAb
KM8014	Anti S1-Tag Mouse mAb
KM8016	Anti Tag-100-Tag Mouse mAb
KM8017	Anti mCherry-Tag Mouse mAb
KM8018	Anti HPC4-Tag Mouse mAb
KM8019	Anti SUMO-Tag Mouse mAb
KM8020	Anti RFP-Tag Mouse mAb
LK8001	Anti His-Tag-HRP Mouse mAb
LK8002	Anti Flag-Tag-HRP Mouse mAb
LK8003	Anti c-Myc-Tag-HRP Mouse mAb
LK8009	Anti GFP-Tag-HRP Mouse mAb
LK2003	Goat anti-Mouse IgG(H+L)-HRP
LK2004	Goat anti-Mouse IgG(H+L)-AP
GM200G-02	FITC Goat anti-Mouse IgG(H+L)
GM200G-09	PE Goat anti-Mouse IgG(H+L)
GM200G-37	IF555 Goat anti-Mouse IgG(H+L)
GM200G-38	IF488 Goat anti-Mouse IgG(H+L)
GM200G-40	IF647 Goat anti-Mouse IgG(H+L)
GM200G-43	IF594 Goat anti-Mouse IgG(H+L)

For Research Use Only.