

APC Anti-Mouse CD45R/B220 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

| Catalog Number | Vial Size |
|----------------|-----------|
| M10451-11A | 25 µg |
| M10451-11C | 100 µg |

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

| Clone | Isotype | Reactivity |
|-------------|---------|------------|
| RA3.3A1/6.1 | Rat IgM | Mouse |

Description

CD45R, also known as B220, is an isoform of CD45. It is a member of the protein tyrosine phosphatase (PTP) family with a molecular weight approximately 180-240 kD. CD45R is expressed on B cells (at all developmental stages from pro-B cells through mature B cells), activated B cells, subsets of T and NK cells. CD45R (B220) is also expressed on a subset of abnormal T cells involved in the pathogenesis of systemic autoimmunity in MRL-Fas^{lpr} and MRL-Fas^{gld} mice. It plays a critical role in TCR and BCR signaling. The primary ligands for CD45 are galectin-1, CD2, CD3, and CD4. CD45R is commonly used as a pan-B cell marker; however, CD19 may be more appropriate for B cell specificity.

Illustration of Immunofluorescent Staining

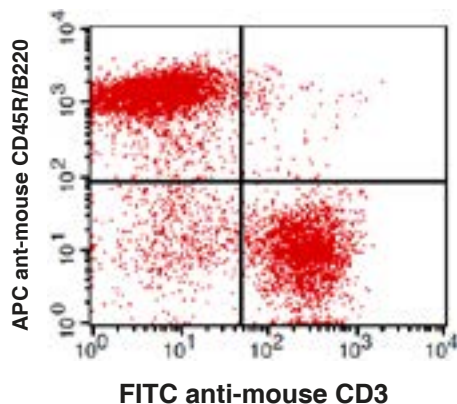

C57BL/6 mouse splenocytes stained with FITC anti-mouse CD3 and APC anti-mouse CD45R/B220

Product Information

Conjugation: APC

Formulation: PBS pH 7.2, 0.09% NaN₃, 0.2% BSA

Concentration: 0.2 mg/ml

Storage: Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used ≤ 0.25 µg /10⁶ cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

References

- [1] Shih, F.F., et al. 2006. J. Immunol. 176:3438.
- [2] Bouwer, H.G.A., et al. 2006. P. Natl. Acad. Sci. USA 103:5102.

For Research Use Only.