

APC Anti-Mouse CD80 (B7-1) Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10801-11A	25 µg
M10801-11C	100 µg

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
16-10A1	Hamster IgG	Mouse

Description

CD80 is a 60 kD highly glycosylated protein. It is a member of the Ig superfamily and is also known as B7-1, B7, and Ly-53. CD80 is constitutively expressed on dendritic cells and monocytes/macrophages, and inducibly expressed on activated B and T cells. The ligation of CD28 on T cells with CD80 and CD86 (B7-2) on antigen presenting cells (such as dendritic cells, macrophages, and B cells) elicits co-stimulation of T cells resulting in enhanced cell activation, proliferation, and cytokine production. CD80 appears to be expressed later in the immune response than CD86. CD80 can also bind to CD152, also known as CTLA-4, to deliver an inhibitory signal to T cells.

Illustration of Immunofluorescent Staining

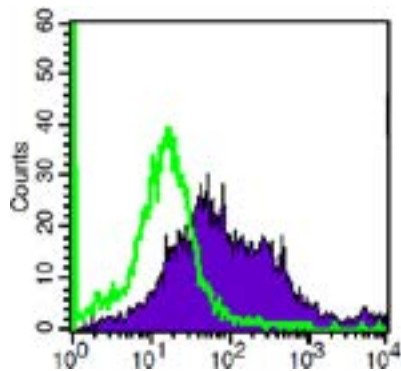

Log Fluorescence Intensity

LPS (3-day) stimulated C57BL/6 mouse splenocytes stained with APC anti-mouse CD80

Product Information

Conjugation: APC

Formulation: PBS pH 7.2, 0.09% NaN₃, 0.2% BSA

Concentration: 0.2 mg/ml

Storage: Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used $\leq 1.0 \mu\text{g} / 10^6$ cells in 100 μl). Since applications vary, the appropriate dilutions must be determined for individual use.

References

- [1] Barclay AN, et al. 1997. The Leukocyte Antigen FactsBook Academic Press.
- [2] Linsley PS, et al. 1991. J. Exp. Med. 174:561.
- [3] Salomon B, et al. 2001. Annu. Rev. Immunol. 19:225.

For Research Use Only.