

## APC Anti-Mouse CD90.2 (Thy-1.2) Monoclonal Antibody



天津三箭生物技术股份有限公司  
Tianjin Sungene Biotech Co., Ltd.  
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10902-11A	25 µg
M10902-11C	100 µg

**Market** | 400-621-0003  
marketing@sungenebiotech.com

**Support** | 022-66211636-8024  
techsupport@sungenebiotech.com

**Web** | www.sungenebiotech.com

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.  
This product is guaranteed up to one year from purchase.


### Purified Antibody Characterization

Clone	Isotype	Reactivity
30H12	Rat IgG2b	Mouse

### Description

CD90.2 is a 25-35 kD immunoglobulin superfamily member also known as Thy1.2. It is expressed on hematopoietic stem cells and neurons, all thymocytes, and peripheral T cells in Thy1.2 bearing mouse strains (Balb/c, CBA/J, C3H/He, C57BL/-, DBA, NZB/-). CD90.2 is a glycosylphosphatidylinositol (GPI)-anchored membrane glycoprotein involved in signal transduction. CD90.2 is involved in costimulation of lymphocyte proliferation and induction of hematopoietic stem cells differentiation. CD90.2 has been shown to interact with CD45. The 30H12 antibody has been reported to induce Ca<sup>2+</sup> flux in thymocytes and, in combination with antibody against the CD3/TCR complex, promote thymocyte apoptosis and inhibit CD3-mediated proliferative responses of mature T lymphocytes.

### Illustration of Immunofluorescent Staining

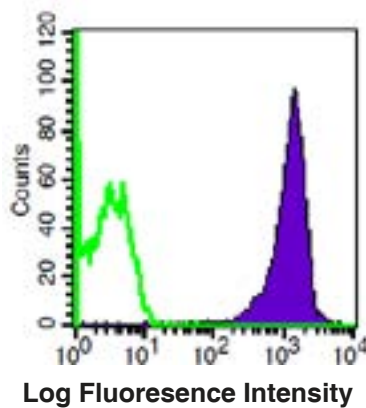

C57BL/6 mouse thymocytes stained with APC anti-mouse CD90.2

### Product Information

**Conjugation:** APC

**Formulation:** PBS pH 7.2, 0.09% NaN<sub>3</sub>, 0.2% BSA

**Concentration:** 0.2 mg/ml

**Storage:** Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

**Application:** Recommended Application: FC

**Usage:** Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used ≤ 0.25 µg /10<sup>6</sup> cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

### References

- [1] Barclay A, et al. 1997. The Leukocyte Antigen FactsBook Academic Press.
- [2] Craig W, et al. 1993. J. Exp. Med. 177:1331.
- [3] Reif AE and Schlesinger M. 1989. Cell Surface Antigen Thy-1.
- [4] Mayani H, et al. 1994. Blood 83:2410.

**For Research Use Only.**