

Biotin Anti-Mouse CD119 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M11192-08B	50 µg
M11192-08E	500 µg

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
GR-20	Rat IgG2a	Mouse

Description

CDw119 is a 90 kD immunoglobulin superfamily member, also known as IFN- γ R α chain. It is a class II cytokine receptor family member that serves as a IFN- γ -binding chain associated with the IFN- γ β chain also known as AF-1. In addition to ligand binding, CDw119 participates in ligand trafficking. CDw119 is expressed on T and B cells, NK cells, fibroblasts, endothelial, and epithelial cells. Binding of IFN- γ induces receptor dimerization, internalization, Jak1 and Jak2 protein kinase activation and, ultimately, STAT1 activation. IFN- γ initiates and regulates a variety of immune responses.

Illustration of Immunofluorescent Staining

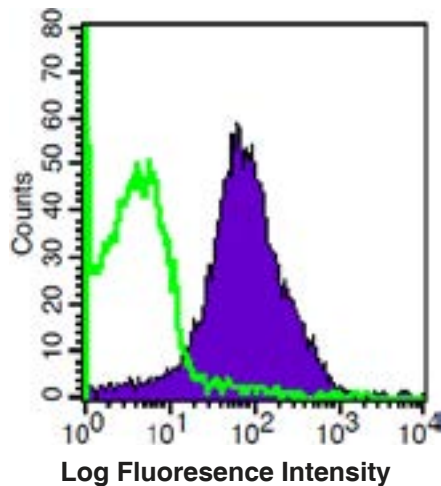

C57BL/6 mouse splenocytes stained with Biotin anti-mouse CD119, followed by PE-SA

Product Information

Conjugation: Biotin

Formulation: PBS pH 7.2, 0.09% NaN₃, 0.2% BSA

Concentration: 0.5 mg/ml

Storage: Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used ≤ 0.25 µg /10⁶ cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

References

[1] Bach EA, et al. 1995. Science 270: 1215.

For Research Use Only.