

## FITC Anti-Mouse CD40 Monoclonal Antibody



天津三箭生物技术股份有限公司  
Tianjin Sungene Biotech Co., Ltd.  
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10401-02B	50 µg
M10401-02E	500 µg

**Market** | 400-621-0003  
marketing@sungenebiotech.com

**Support** | 022-66211636-8024  
techsupport@sungenebiotech.com

**Web** | www.sungenebiotech.com

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.  
This product is guaranteed up to one year from purchase.


### Purified Antibody Characterization

Clone	Isotype	Reactivity
FGK4.5	Rat IgG2a	Mouse

### Description

CD40 is a 48 kD type I transmembrane glycoprotein also known as Bp50. It is a member of the tumor necrosis factor receptor (TNFR) superfamily and is expressed on B cells, basal epithelial cells, macrophages, follicular dendritic cells, endothelial cells, and a subset of CD34+ hematopoietic progenitors. CD40 regulates B cell development/maturation, Ig isotype switching and, in combination with other signals such as IL-4, protects B cells from surface Ig-induced apoptosis and promotes proliferation. Interaction of CD40 with its ligand CD154(gp39), which is expressed on activated T cells, is important in costimulation and immune regulation.

### Illustration of Immunofluorescent Staining

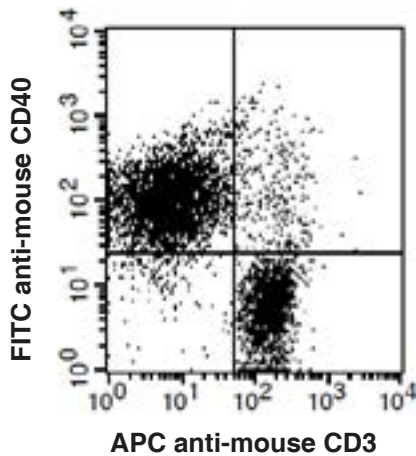

C57BL/6 mouse splenocytes stained with FITC anti-mouse CD40 and APC anti-mouse CD3

### Product Information

**Conjugation:** FITC

**Formulation:** PBS pH 7.2, 0.09% NaN<sub>3</sub>, 0.2% BSA

**Concentration:** 0.5 mg/ml

**Storage:** Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

**Application:** Recommended Application: FC

**Usage:** Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used  $\leq 1.0 \mu\text{g} / 10^6$  cells in 100  $\mu\text{l}$ ). Since applications vary, the appropriate dilutions must be determined for individual use.

### References

- [1] Barclay A, et al. 1997. The Leukocyte Antigen FactsBook Academic Press.
- [2] Bancherou J, et al. 1994. Annu. Rev. Immunol. 12:881.
- [3] Clark EA, et al. 1996. P. Natl. Acad. Sci. USA 83:4494.

**For Research Use Only.**