

IF555 Goat Anti-Rabbit IgG(H+L)

天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Catalog Number	Vial Size
GR200G-37C	100 ug
GR200G-37E	500 ug

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.

Purified Antibody Characterization

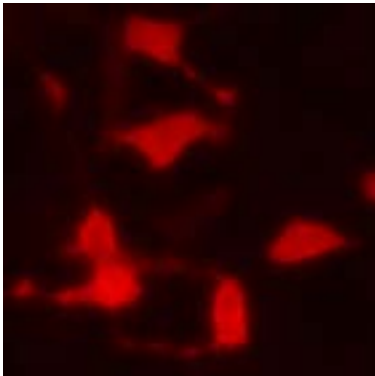
Isotype	Reactivity
Polyclonal Goat IgG	Rabbit

Description

Preparation: The polyclonal antibody was purified from goat antiserum by rabbit IgG affinity chromatography.

Purity: >90% by SDS-PAGE

Specificity: The antibody reacts with whole molecule Rabbit IgG. No antibody was detected against non-immunoglobulin serum proteins. The antibody may cross-react with immunoglobulins from other species.

Illustration of Immunofluorescent Staining

Immunofluorescence Staining of HEK293 cell transfected with Myc-tag fusion plasmid. The primary antibody is rabbit polyclonal antibody. The secondary antibody is IF555 goat anti-rabbit IgG(H+L) diluted at 1/200.

Product Information

Conjugation: IF555

Formulation: PBS pH 7.2, 0.09% NaN₃, 0.2% BSA

Concentration: 1.0 mg/ml

Storage: Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: IF(Immunofluorescent). Not yet tested in other applications.

Usage: IF: 1/100 ~ 1/200

Since applications vary, the appropriate dilutions must be determined for individual use.

iFluor™555 (IF555) has fluorescence excitation and emission maxima close to 555 nm and 570 nm. These spectral characteristics make it excellent alternatives to Cy3® and Alexa Fluor® 555.

For Research Use Only.