

## PE Anti-Human CD22 Monoclonal Antibody



天津三箭生物技术股份有限公司  
Tianjin Sungene Biotech Co., Ltd.  
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
H20221-09G	25 tests
H20221-09H	100 tests

<b>Market</b>	400-621-0003 marketing@sungenebiotech.com
<b>Support</b>	022-66211636-8024 techsupport@sungenebiotech.com
<b>Web</b>	www.sungenebiotech.com

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.  
This product is guaranteed up to one year from purchase.


### Purified Antibody Characterization

Clone	Isotype	Reactivity
HIB22	Mouse IgG1	Human

### Description

CD22 is a 130 kD type I transmembrane glycoprotein also known as Siglec-2 and BL-CAM. It is a member of the immunoglobulin superfamily (sialoadhesion subgroup). CD22 is expressed in the cytoplasm of pro-B and pre-B cells, and on the surface of mature B and activated B cells, but not on plasma cells. CD22 is present in the B cell receptor complex and associates with SHP-1, Syk, Lck, Lyn, and phospholipase C $\gamma$ 1. A primary function of CD22 is thought to be in limiting antigen receptor signaling by modulating B cell activation threshold. CD22 has been shown to bind to CD45RO and CD75, although the natural ligands for this molecule remain controversial

### Illustration of Immunofluorescent Staining

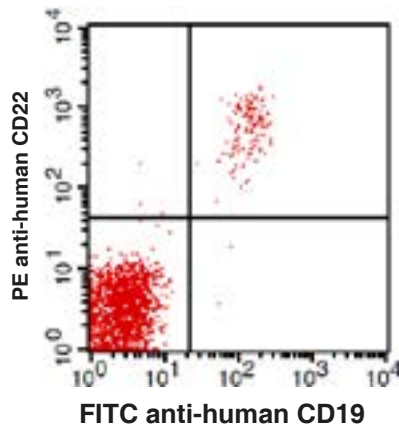

Human peripheral blood lymphocytes stained with FITC anti-human CD19 and PE anti-human CD22

### Product Information

**Conjugation:** PE

**Formulation:** PBS pH 7.2, 0.09% NaN<sub>3</sub>, 0.2% BSA

**Storage:** Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

**Application:** Recommended Application: FC

**Usage:** Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used from 20  $\mu$ L to 5  $\mu$ L per 100  $\mu$ L of peripheral blood. Please check your vial). Since applications vary, the appropriate dilutions must be determined for individual use.

### References

- [1] Schlossman, S., et al. Eds. 1995. Leukocyte Typing V: White Cell Differentiation Antigens. Oxford University Press. New York.
- [2] Clark, E., 1993. J. Immunol.. 150:4715.
- [3] Shan, D. and O. Press. 1995. J. Immunol.. 154:4466.

**For Research Use Only.**