

PE Anti-Human CD29 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
标准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
H20291-09G	25 tests
H20291-09H	100 tests

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
HI29a	Mouse IgG1	Human

Description

CD29 is a 130 kD single chain type I glycoprotein, known as integrin $\beta 1$, VLA- β chain, or gpIIa. It is broadly expressed on a majority of hematopoietic and non-hematopoietic cells, including leukocytes (although at low level on granulocytes), platelets, fibroblasts, endothelial cells, epithelial cells, and mast cells. CD29 is a member of the integrin family. It is non-covalently associated with integrin $\alpha 1$ - $\alpha 6$ chains to form VLA-1 to VLA-6 molecules, respectively. Integrins which include CD29 bind to several cell surface (e.g. VCAM-1, MadCAM-1) and extracellular matrix molecules. CD29 acts as a fibronectin receptor and is involved in a variety of cell-cell and cell-matrix interactions.

Illustration of Immunofluorescent Staining

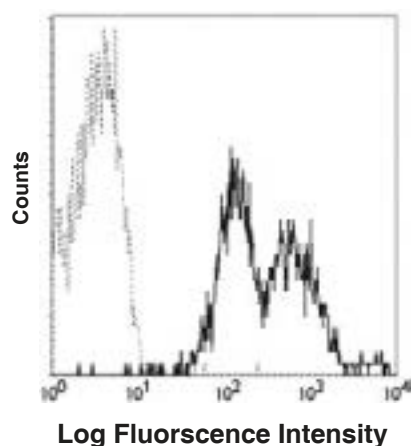

Human peripheral blood lymphocytes
stained with PE anti-human CD29

Product Information

Conjugation: PE

Formulation: PBS pH 7.2, 0.09% NaN_3 , 0.2% BSA

Storage: Keep as concentrated solution.
Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used from 20 μL to 5 μL per 100 μL of peripheral blood. Please check your vial). Since applications vary, the appropriate dilutions must be determined for individual use.

References

- [1] Hemler M. 1990. Annu. Rev. Immunol. 8:365.
- [2] Hynes R. 1992. Cell 69:11.

For Research Use Only.