

PE Anti-Human CD54 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
H10541-09G	25 tests
H10541-09H	100 tests

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
R6-5-D6	Mouse IgG2a	Human

Description

CD54 is a 85-110 kD type I transmembrane protein also known as ICAM-1. It is expressed on activated endothelial cells, high endothelial venules, T and B cells, monocytes/macrophages, granulocytes, and dendritic cells. The expression of ICAM-1 can be released from the cell surface. CD54 plays a role in cellular adhesion and is involved in inflammation and leukocyte extravasation. CD54 has also been shown to be the major cellular receptor for rhinovirus. ICAM-1 binds to CD11a/CD18 (LFA-1), CD11b/CD18 (Mac-1), CD11c/CD18 (p150, 95) as well as hyaluronan and fibrinogen.

Illustration of Immunofluorescent Staining

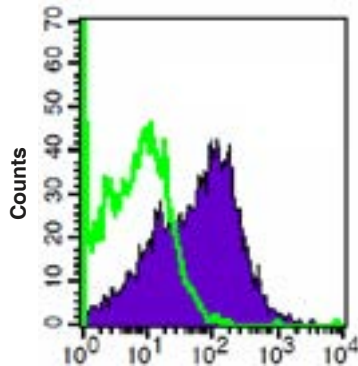

Log Fluorescence Intensity

Human peripheral blood lymphocytes stained
with PE anti-human CD54

Product Information

Conjugation: PE

Formulation: PBS pH 7.2, 0.09% NaN₃,
0.2% BSA

Storage: Keep as concentrated solution.
Store at 4°C and protected from prolonged
exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality
control tested by immunofluorescent staining
with flow cytometric analysis (The amount
of the reagent is suggested to be used from
20 μL to 5 μL per 100 μL of peripheral blood.
Please check your vial). Since applications
vary, the appropriate dilutions must be
determined for individual use.

References

- [1] Voraberger G, et al. 1991. J. Immunol. 147:2777.
- [2] Staunton DE, et al. 1988. Cell 52:925.
- [3] Greve JM, et al. 1989. Cell 56:839.

For Research Use Only.