

PE/Cy5 Anti-Mouse CD16/32 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10161-35A	25 µg
M10161-35C	100 µg

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
2.4G2	Rat IgG2b	Mouse

Description

The 2.4G2 antibody reacts specifically with a common nonpolymorphic epitope on the extracellular domains of the mouse FcγIII and FcγII receptors. It has also been reported to bind the FcγI receptor (CD64) via its Fc domain. 2.4G2 mAb blocks non-antigen-specific binding of immunoglobulins to the FcγIII and FcγII, and possibly FcγI, receptors in vitro and in vivo. CD16 and/or CD32 are expressed on natural killer cells, monocytes, macrophages, dendritic cells (at low levels), Kupffer cells, granulocytes, mast cells, B lymphocytes, immature thymocytes, and some activated mature T lymphocytes.

Illustration of Immunofluorescent Staining

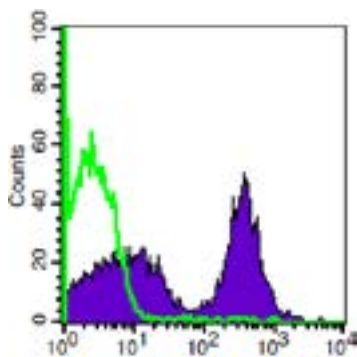

Log Fluorescence Intensity

C57BL/6 mouse splenocytes stained with
PE/Cy5 anti-mouse CD16/32

Product Information

Conjugation: PE/Cy5

Formulation: PBS pH 7.2, 0.09% NaN₃,
0.2% BSA

Concentration: 0.2 mg/ml

Storage: Keep as concentrated solution.
Store at 4°C and protected from prolonged
exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used ≤ 0.25 µg /10⁶ cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

References

- [1] Araujo-Jorge T, Rivera MT, el Bouhdidi A, Daeron M, Carlier Y. An Fc gamma RII-, Fc gamma RIII-specific monoclonal antibody (2.4G2) decreases acute Trypanosoma cruzi infection in mice. *Infect Immun.* 1993; 61(11):4925-4928.(Clone-specific)
- [2] Benhamou M, Bonnerot C, Fridman WH, Daeron M. Molecular heterogeneity of murine mast cell Fc gamma receptors. *J Immunol.* 1990; 144(8):3071-3077.
- [3] Mellman IS, Unkeless JC. Purification of a functional mouse Fc receptor through the use of a monoclonal antibody. *J Exp Med.* 1980; 152(4):1048-1069.

For Research Use Only.