

PE/Cy5 Anti-Mouse CD28 Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10282-35A	25 µg
M10282-35C	100 µg

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
37.51	Hamster IgG	Mouse

Description

CD28 is a 44 kD glycoprotein, also known as Tp44 or T44. It is a member of the Ig superfamily, expressed on thymocytes, most peripheral T cells, and NK cells. In association with CD80 (B7-1) and CD86 (B7-2), CD28 acts as the second signal for T and NK cell activation and proliferation. The 37.51 antibody has been reported to augment in vitro T cell proliferation and cytokine production, and promote CTL development.

Illustration of Immunofluorescent Staining

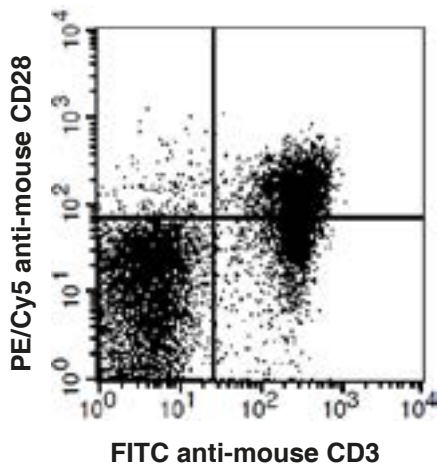

C57BL/6 mouse splenocytes stained with
FITC anti-mouse CD3 and PE/Cy5 anti-mouse CD28

Product Information

Conjugation: PE/Cy5

Formulation: PBS pH 7.2, 0.09% NaN₃,
0.2% BSA

Concentration: 0.2 mg/ml

Storage: Keep as concentrated solution.
Store at 4°C and protected from prolonged
exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality
control tested by immunofluorescent staining
with flow cytometric analysis (The amount of
the reagent is suggested to be used ≤ 1.0
µg /10⁶ cells in 100 µl). Since applications
vary, the appropriate dilutions must be
determined for individual use.

References

- [1] Barclay AN, et al. 1997. The Leukocyte Antigen FactsBook Academic Press.
- [2] Lenschow DJ, et al. 1996. Annu. Rev. Immunol. 14:233.
- [3] Gross JA, et al. 1992. J. Immunol. 149:380.

For Research Use Only.