

# PerCP Anti-Rat CD4 (domain 1) Monoclonal Antibody



天津三箭生物技术股份有限公司  
Tianjin Sungene Biotech Co., Ltd.  
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
R10041-32A	25 µg
R10041-32C	100 µg

**Market** | 400-621-0003  
marketing@sungenebiotech.com

**Support** | 022-66211636-8024  
techsupport@sungenebiotech.com

**Web** | www.sungenebiotech.com

**Important Note:** Centrifuge before opening to ensure complete recovery of vial contents.  
This product is guaranteed up to one year from purchase.


## Purified Antibody Characterization

Clone	Isotype	Reactivity
OX-38	Mouse IgG2a	Rat

## Description

CD4, also known as T4, is a 55kD glycoprotein member of the immunoglobulin superfamily and is expressed on majority of thymocytes, macrophages, and a peripheral T cell subset (T helper cells). CD4 is a T cell co-receptor that interacts with the MHC class II molecule and is involved in T cell activation.

## Illustration of Immunofluorescent Staining

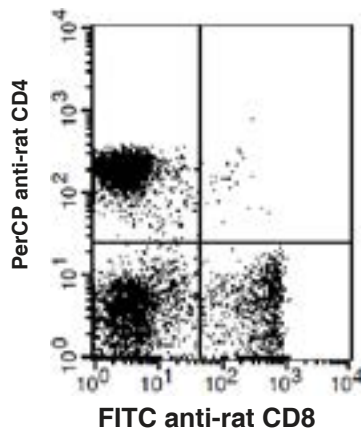

Wister rat splenocytes stained with PerCP anti-rat CD4 and FITC anti-rat CD8 antibody

## Product Information

**Conjugation:** PerCP

**Formulation:** PBS pH 7.2, 0.09% NaN<sub>3</sub>, 0.2% BSA

**Concentration:** 0.2 mg/ml

**Storage:** Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

**Application:** Recommended Application: FC

**Usage:** Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used ≤ 0.25 µg /10<sup>6</sup> cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

## References

[1] Pettigrew, G.J., et al. 1998. J. Immunol. 161:1292.

**For Research Use Only.**