

PerCP/Cy5.5 Anti-Mouse CD4 (YTS191) Monoclonal Antibody



天津三箭生物技术股份有限公司
Tianjin Sungene Biotech Co., Ltd.
精准 高效 稳定 Precision Efficient Stable

Catalog Number	Vial Size
M10041-33A	25 µg
M10041-33C	100 µg

Market | 400-621-0003
marketing@sungenebiotech.com

Support | 022-66211636-8024
techsupport@sungenebiotech.com

Web | www.sungenebiotech.com

Important Note: Centrifuge before opening to ensure complete recovery of vial contents.
This product is guaranteed up to one year from purchase.


Purified Antibody Characterization

Clone	Isotype	Reactivity
YTS191	Rat IgG2b	Mouse

Description

CD4 is a 55 kD protein also known as L3T4 or T4. It is a member of the Ig superfamily, primarily expressed on most thymocytes, a subset of T cells, and weakly on macrophages and dendritic cells. It acts as a coreceptor with the TCR during T cell activation and thymic differentiation by binding MHC class II and associating with the protein tyrosin kinase, lck.

Illustration of Immunofluorescent Staining

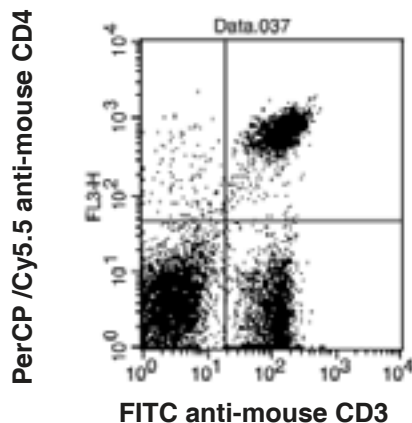

C57BL/6 mouse splenocytes stained with FITC anti-mouse CD3 and PerCP/Cy5.5 anti-mouse CD4

Product Information

Conjugation: PerCP/Cy5.5

Formulation: PBS pH 7.2, 0.09% NaN₃, 0.2% BSA

Concentration: 0.2 mg/ml

Storage: Keep as concentrated solution. Store at 4°C and protected from prolonged exposure to light. **Do not freeze.**

Application: Recommended Application: FC

Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis (The amount of the reagent is suggested to be used $\leq 0.25 \mu\text{g} / 10^6$ cells in 100 µl). Since applications vary, the appropriate dilutions must be determined for individual use.

References

- [1] Barclay, A., et al. 1997. The Leukocyte Antigen FactsBook, Academic Press.
- [2] Bierer, B.E., et al. 1989. Annu. Rev. Immunol. 7:579.
- [3] Janeway, C.A. 1992. Annu. Rev. Immunol. 10:645.

For Research Use Only.